


SOLID STATE HRM CATHETERS

A FAST AND RELIABLE WAY FOR YOUR MANOMETRY TESTS

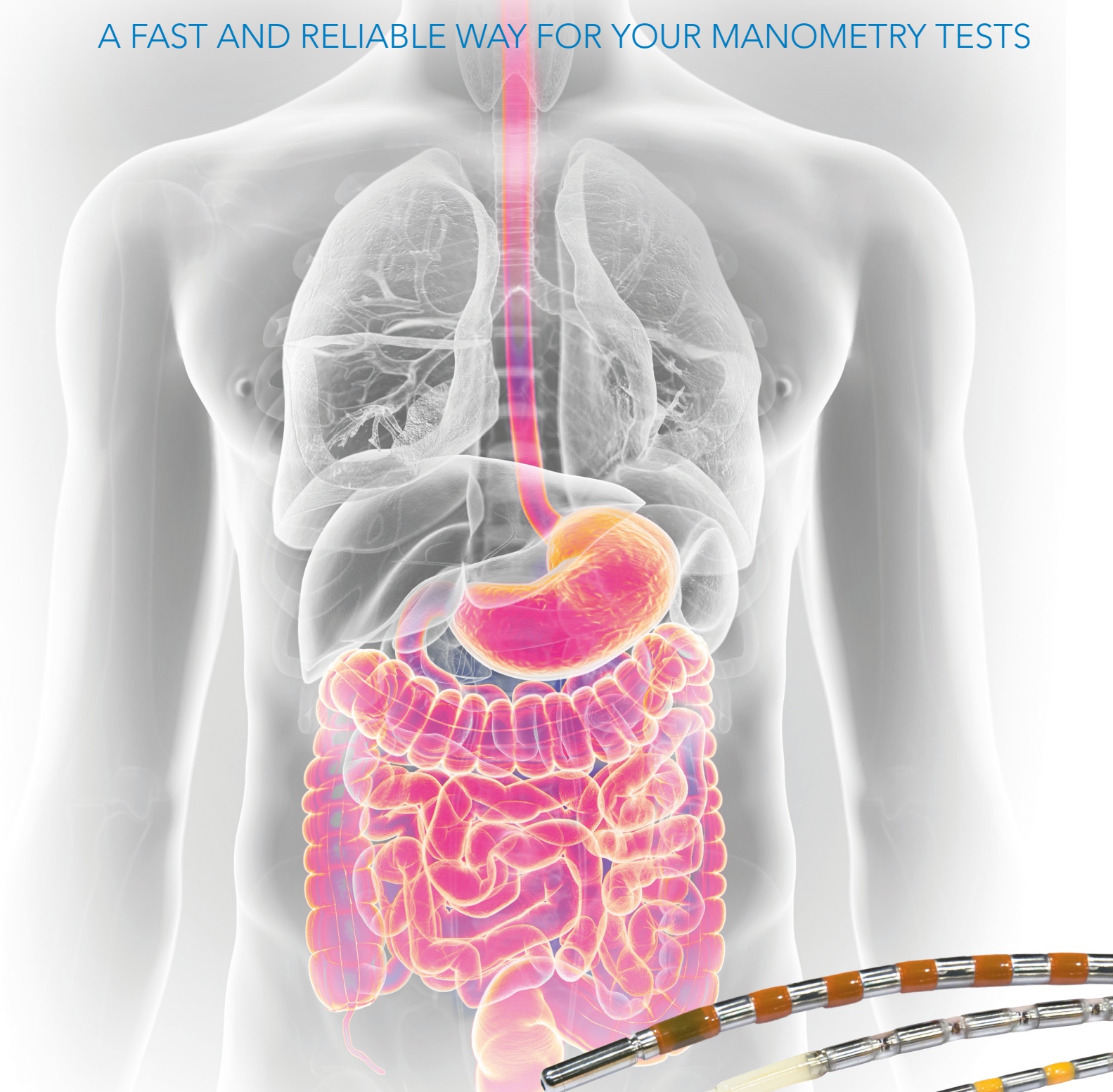

CHOICE
EASE OF USE
RELIABILITY


SOLID STATE HRM CATHETERS

Solid state manometry has been a proven technique for over 15 years that offers rapid response to measure the true physiological signature. Reusable solid state catheters need no study preparation, no calibration and use only one connector for all pressures.

LABORIE is offering a full range of solid state catheters for:

- HRM and HRIM (esophageal)
- HRAM and 360 degrees HRAM (anorectal)
- HRCM (colon)
- HRSM (small bowel/antroduodenum)
- Sphincter of Oddi

in a variety of diameters ranging from 6 French up to 14 French, to ensure accurate measurement for different patient groups. Our catheters are ranging from 1 pressure for simple ARM studies up to 36 pressure channels for HRM and 40 pressure channels for 360 degrees HRAM. Optionally, impedance electrodes and lumen can be added.

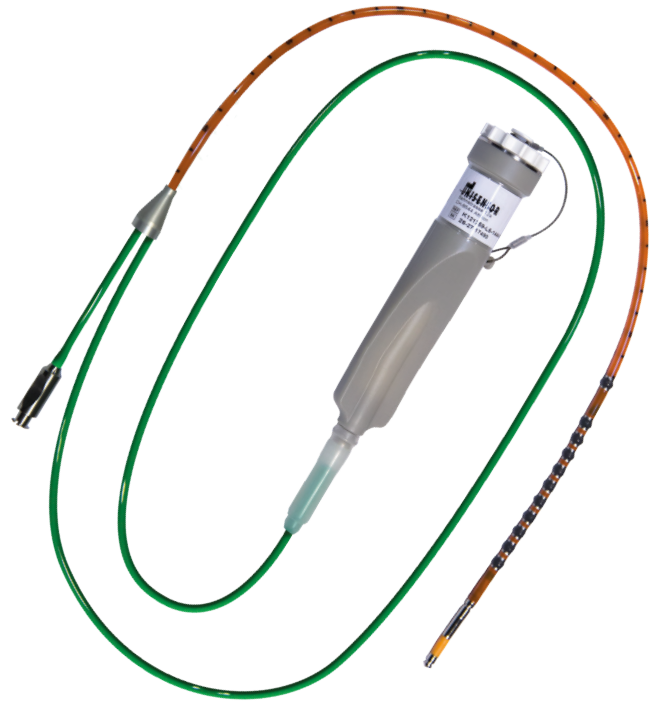

WHAT DOES THIS MEAN FOR YOU AND YOUR PATIENTS?

CHOICE


You have a great variety of standard and special-order catheters to choose from, offering a solution for each of your patient groups, each of your investigation types and for each budget. We could even make a customized catheter to your own design for research or for very specific investigations, on request.

EASE OF USE

Our solid state catheters need no study preparation or calibration and use only one connector for all pressures. The catheters can be cleaned mechanically, rather than manually, due to the special watertight cap. All this allows you to spend more time with your patients.

RELIABILITY

The high quality and the minimal drift of our catheters make them very suitable for your prolonged investigations, such as colonic and AD manometry studies. Our Premium service contract will allow you to perform your investigations without any interruptions, even in case of a repair, and will help you to get the most out of your catheter.

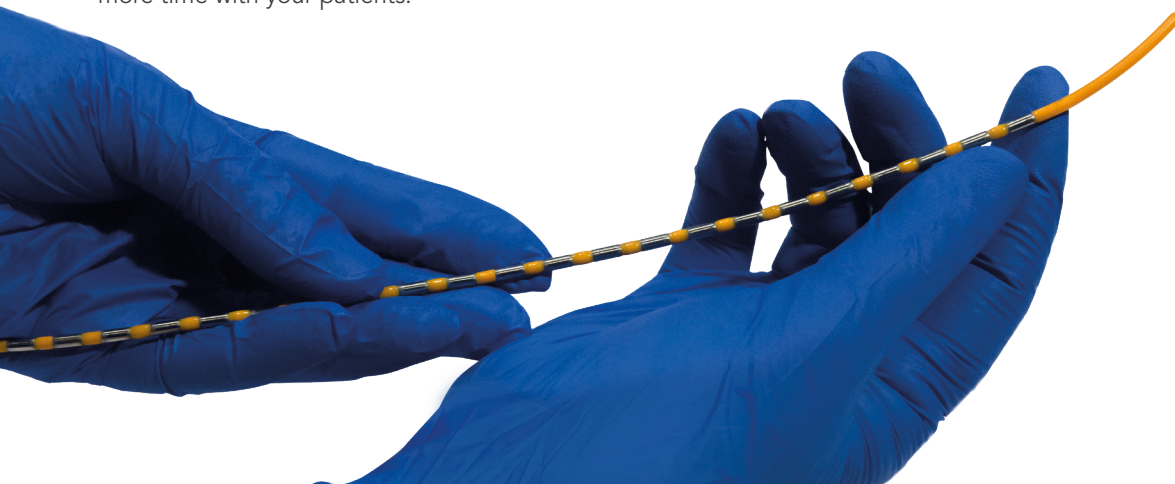

STANDARD SOLID STATE HRM CATHETERS

Sales code	Investigation	# pressures	Type of pressures	Diameter	Impedance
K12959-L5-1038-D	HRAM	9	1 for balloon, 8 circumferential	12 French	-
K121259-L5-1444-D	HRAM	12	1 for balloon, 11 circumferential	12 French	-
K122359-L5-1623-D	HRAM (3D)	23	All unidirectional, 5x4 at 90 degrees	12 French	-
K102459-00-1389-D	HRM	24	All unidirectional, P2 at 180 degrees	10 French	-
K103659-00-1102-D	HRM	36	5 circular (LES), 31 unidirectional	10 French	-
K83659-E-10107-D	HRIM	36	All unidirectional	8 French	12Z (13 rings)
K103659-E-1545-D	HRIM	36	5 circular (LES), 31 unidirectional	10 French	16Z (17 rings)

Standard catheters: < 2 weeks delivery time, catheter definitions standard included in the software.

SPECIAL-ORDER SOLID STATE HRM CATHETERS

Sales code	Investigation	# pressures	Type of pressures	Diameter	Impedance
K123659-00-1002-D	HRM	36	All circumferential	12 French	-
K62559-00-1049-D	HRM	25	All unidirectional	6 French	-
K62559-E-1022-D	HRIM	25	All unidirectional	6 French	12Z (13 rings)
K143659-L6-1375-D	HRM Colonic	36	All unidirectional	14 French	-
K143659-L6-1352-D	HRM Colonic	36	All unidirectional	14 French	-
K143659-L6-1612-D	HRM AD	36	All unidirectional	14 French	-
K103659-00-1744-D	HRM AD	36	All unidirectional	10 French	-

Special-order catheters: delivery times 6-8 weeks; catheter definitions standard included in the software.
Customized designs are possible on request.



SOLAR GI HRM

THE MOST ADVANCED HRM SYSTEM ON THE MARKET

Our solid state catheters are best used with the Solar GI HRM system, currently the most advanced HRM system on the market, that gives you the ability to:

- Perform quicker (10 minutes only), more reliable and more accurate procedures
- Do better, easier and quicker data analysis
- Lower your procedure costs

PREMIUM SERVICE CONTRACT FOR SOLID STATE HRM CATHETERS ⁽¹⁾

Our optional Premium service contract for solid state HRM catheters includes:


- Extended coverage and warranty
- A replacement loaner within 48-72 hours
- Reduced repair costs

Type of service	Standard warranty	Premium service contract
Number of procedures	200	300
Check-up of the catheter and recalibration	Not included	Included
Loaner replacement	No	48-72 hours
Repairs	Not included	Included
Professional cleaning	Not included	Included
Accidental damage	Not included	50% off cost

The premium service contract will allow you to perform your motility tests with no disruptions in your practice for the duration of repairs. At the same time you will benefit from an increased number of uses under warranty for your catheter.

The Premium service contract is available for newly purchased STANDARD catheters only.

⁽¹⁾ Our Premium service contract may not be available in certain regions. Contact your sales representative for local availability.

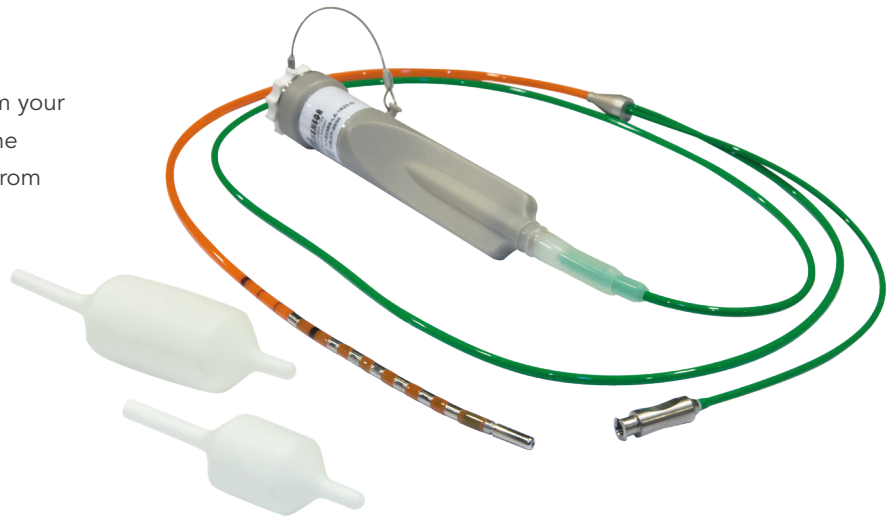

Data can be subject to change without notice 11-2019 © Copyright by Laborie The Netherlands 0077-BRG-016-EN-LTR-V1.01 Solid State HRM catheters PRC-023680



GI OFFICE:
Tel: +31 20 893 80 80
Email: GImarketing@laborie.com

De Cuserstraat 93
1081 CN Amsterdam
The Netherlands

USA:
Tel.: +1 800 522 6743
Email: usmarketing@laborie.com

www.laborie.com